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INSTRUCTIONS FOR THE USE OF PARALLELING PINS

1. Product description

Reusable invasive medical grade titanium surgical device for temporary use as a guide to determine drill depth levels and angles before proceeding with the drilling for a second implant during the same surgical procedure.

2. Planned use:

Paralleling pins are used to determine the parallelism between teeth and the structures supporting the final prosthesis.

3. Directions for use

Paralleling pins are used to determine the parallelism between teeth and the structures supporting the final prosthesis.

4. Contraindications

Allergy and hypersensitivity to the chemical components of the materials used (medical grade titanium).

5. Side effects and interactions

Follow the instructions for use recommended by the manufacturer. Incorrect use of the product may cause injury to the patient.

6. Cleaning, disinfection and sterilisation

In the event of any packaging defect, the product should not be used and the manufacturer notified immediately.

All BIONER paralleling pins are supplied clean and appropriately packaged. They require sterilisation by the healthcare professional before use in the mouth.

Bioner SA always recommends performing previous cleaning with products commonly used in dentistry, with emphasis on ultrasonic baths or automatic cleaning and disinfection methods before proceeding with sterilisation.

Remove the product from its original packaging and place it in a suitable container before sterilisation.

Sterilisation should be performed in sterilisation bags or boxes in a wet steam autoclave at 121°C for 15 minutes or 134°C for at least 3 minutes, as per UNE-EN-ISO 17665-1:2007.

Do not use the product after sterilisation if there is any detectable deterioration of the container.

7. Warnings

Paralleling pins may be reused after cleaning, disinfection and sterilisation.

The healthcare professional using the product must have experience in using this product safely and in accordance with instructions for use provided by BIONER SA, the said professional being liable for the suitability of the component to the patient's needs.

8. Precautions

It is especially important to ensure that each reference is used in combination with dimensionally compatible elements as indicated on the product label.

9. Compatibility


BIONER paralleling pins family references must be used in conjunction with compatible components. Carefully examine the components available in the company product catalogue for suitability before insertion.

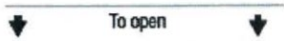10. Validity

These instructions for use correspond to the current version of the paralleling pins family references indicated at the foot of the document.

11. Storage

BIONER products do not require any special storage conditions even though it is recommended not to expose them to conditions that could affect the product, its packaging or label, such as extreme temperature and humidity, knocks, etc.

ISO 15223 Symbol	Description
	Indicates the direction for opening the blister
	BIONER SA logotype
REF	Catalogue reference
LOT AOFAAXXXX	Batch code
	Date of Manufacture
	Refer to the instructions for use
	Do not use if package is damaged
Family	Name of the product family
Mat = Titanium	General indication of the material used
q=1	Indicates the number contained in each blister

	Non-sterile
	QR code (reference and batch)
	Consult warnings and precautions in instructions for use
	Manufacturer
	EC Marking

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